



Mollis Trolley Mattress

For comfort in emergencies



Järven Health Care

Bromsvägen 3 | SE-891 60 Örnsköldsvik | SWEDEN
Telephone: +46 (0)660-22 28 50 | Telefax: +46 (0)660 - 824 10
www.jarven.se | healthcare@jarven.se

Mollis Trolley Mattress (antidecubitus 2 layers)

As the patient often has to remain on the trolley mattress for a long time, it is of great importance he or she is comfortable and that the amount of time does not result in additional complications.

The trolley mattress provides pressure relief and can hence reduce the risk of the patient developing pressure ulcers. Mattress is equipped with a durable surface layer that is comfortable for the patient to lie on.

The underside surface slides easily on most types of surfaces. Along the sides of the mattress are the standard six handles that facilitate lateral movement. We also offer an alternate model with three fabric handles on each side of the mattress.

The trolley mattress is symmetrically shaped.

Foam core

Due to the composition of two different foam types, mattress is pressure relieving and comfortable to lie on even for a long time. The flex foam at the top of the mattress is smooth and soft. The mattress is 10 cm high and consists of two layers, 5 cm highly elastic cold foam at the bottom and 5 cm soft flex foam at the top.

Alternative materials

The Mollis trolley mattress is also available with a welded PVC cover. This material somewhat reduces the comfort properties and do not achieve the same high level of comfort.

Mollis Direct with a firmer mattress core is available as an option. Please contact us for further information.



Please read the instruction manual before using the product.

| Item no. | Description | Size |
|----------|---------------------|------------------|
| 2380 | PU 3 fabric handles | 200 x 60 x 10 cm |
| 2380-70 | PU 3 fabric handles | 200 x 70 x 10 cm |
| 2380-PVC | PU 3 fabric handles | 200 x 60 x 10 cm |
| 2382 | PU 6 welded handles | 200 x 60 x 10 cm |
| 2380-70 | PU 6 welded handles | 200 x 70 x 10 cm |
| 2382-PVC | PU 6 welded handles | 200 x 60 x 10 cm |
| 23813 | Flex pillow | 47 x 26 x 9 cm |
| 23810 | Standard pillow | 50 x 60 cm |

Warranty

2 year for manufacturing defects. NOTE: The warranty does not apply if damage to the mattress or pillow is caused by mechanical impact from pens or other sharp objects.

User weight


0-130 kg.

Fire protection

Meets the requirements of SS 876 00 01

Recycling

The mattress is fully recyclable and should be sorted for energy recovery. No dismantling is required.

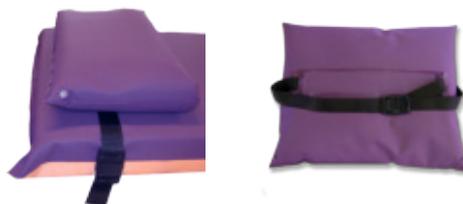

Standard pillow (see picture) and flex pillow (custom-made) are optional.

Results of pressure measurement

| Mattress | Large indenter 300 N | Large indenter 600 N | Small indenter 40 N |
|----------------------------------|----------------------|----------------------|---------------------|
| Mollis Trolley Mattress 2 layers | 36.0 | 73.9 | 67.7 |

A low value means a low pressure against the body, all values in mm Hg. For heavier users it is especially important that the pressure shown under "large indenter, 600 N" is as low as possible.

FACTS ABOUT PRESSURE MEASUREMENT

To bypass the individual impact from different people, Swerea IVF has developed pressure indenters that simulate the shapes of selected parts of the human body. By applying a load equal to that of a human body of average weight on these indenters, it is possible to objectively evaluate different mattresses and beds using pressure mats. This method of measurement is described in the Swedish standard with the designation SS 876 00 13




Järven Health Care

Bromsvägen 3 | SE-891 60 Örnsköldsvik | SWEDEN
 Telephone: +46 (0)660-22 28 50 | Telefax: +46 (0)660 - 824 10
 www.jarven.se | healthcare@jarven.se